

ABLE ASSESS · QUICK START

The mobile platform for **functional assessments.**

Set up Able Assess and run your first standardised falls-risk screening. For the full step-by-step, see the Able Assess Full User Guide v1.10.1.

FOR

Caregivers and assessors
using Able Assess

VERSION

1.10.1


SUPPORT

support@able-care.co

— ONE-TIME SETUP


Get set up — one-time.

Five short steps from a fresh phone to a calibrated, ready-to-use GripAble.

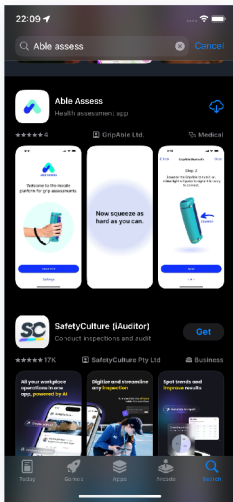
1


Install

Find Able Assess on the App Store or Google Play and download.

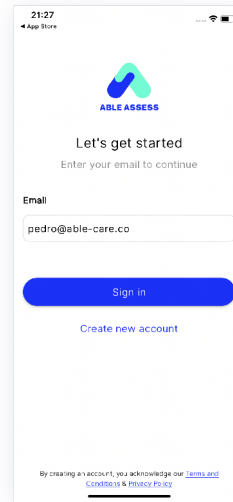
2


Sign in or create account

Open the app. Enter your email, then tap **Sign in** or **Create new account**. You'll receive a one-time-passcode (OTP) by email — paste it in.

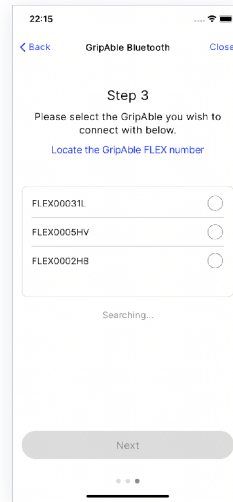
3

Pair your GripAble

Allow Bluetooth, squeeze the GripAble (blue light pulses), then tap its **FLEX number** to pair. The last three characters are unique.

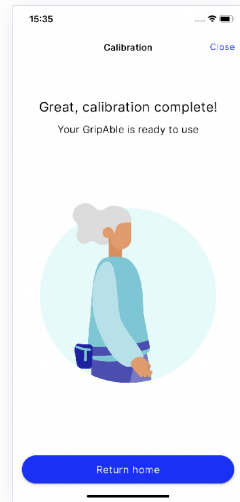
4

Calibrate

Place the GripAble upright on a flat surface and follow the on-screen prompts.

Tip — stay signed in: Enable PIN or biometrics on your phone (HIPAA-compliant). After your first sign-in, just close and reopen the app — you'll skip the OTP step.

Tip — Org account: Email hello@able-care.co if you'd like a shared organisation account set up by our team.

Quick reference

EQUIPMENT

paired and calibrated GripAble · the bag and wrist lanyard · tape measure · sturdy chair without arms · 6m clear walkway

ASSESSMENTS

Single Maximum Grip Strength · Chair Stand Test (CST) · Timed Up and Go (TUG) · 4 Metre Gait Speed (4MGS)

CONTACT


support@able-care.co · we reply within one working day

— EVERYDAY FLOW

Run a Falls Screening — everyday flow.

The full assessment battery. Takes about 5–10 minutes per user.


Before you start: mark a clear walkway with tape at 3m (10ft) for the Timed Up and Go and 6m (20ft) for the 4 Metre Gait Speed.

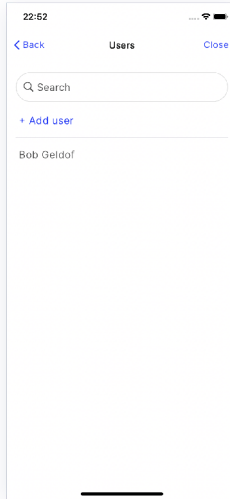

1

Select a user


Tap **Go** to **users**, then tap the user you want to assess (or tap **Add user** to create one).

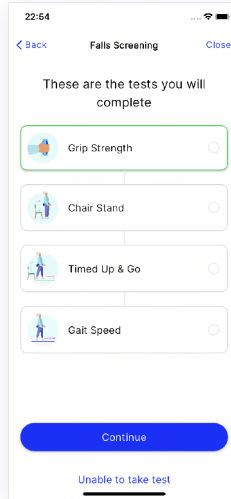

2

Start Falls Screening


From the assessments screen, tap **Falls Screening**.

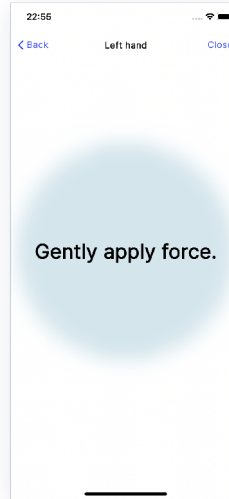

3

Run all four tests


Grip Strength → Chair Stand → Timed Up and Go → 4 Metre Gait Speed. Follow the on-screen setup and instructions for each.

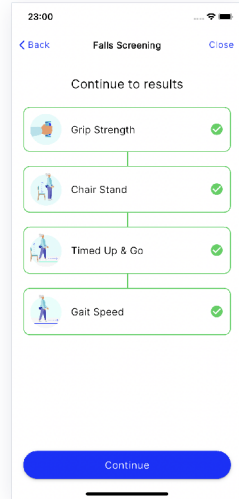

4

Confirm complete


Green ticks across all four tests confirm the screening is complete. Skipped assessments will be greyed out.

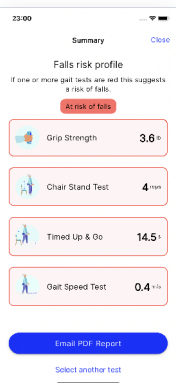

5

Review results

Top-line score per test plus an overall falls-risk indicator. Tap any test to expand and see percentile and threshold.

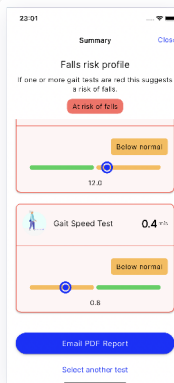

6

Email the PDF report

Tap **Email report** to send a PDF directly to the email of the account you are logged in with. Results are also saved to the Customer Portal at portal.able-care.io.

Note: Only the full Falls Screening generates a PDF. Individual assessments do not.

Able Assess — Quick Start Guide v1.10.1

For the full step-by-step guide, see the **Able Assess Full User Guide v1.10.1**.